

**P91350Hu01**  
**Adipose Differentiation Related Protein (ADRP)**  
**Organism: Homo sapiens (Human)**  
*Instruction manual*

FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

---

5th Edition (Revised in January, 2013)

**[ DESCRIPTION ]**

**Protein Names:** Adipose Differentiation Related Protein

**Synonyms:** ADRP, PLIN2,ADFP

**Species:** Human

**Size:** 100µg

**Source:** *Escherichia coli*-derived

**Subcellular Location:** Membrane; Peripheral membrane protein.

**[ PROPERTIES ]**

**Residues:** Ala5~Tyr232 (Accession # Q99541), with N-terminal His-Tag.

**Grade & Purity:** >95%, 28kDa as determined by SDS-PAGE reducing conditions.

**Formulation:** Supplied as lyophilized form in PBS, pH 7.4, containing 5% sucrose, 0.01% sarcosyl.


**Endotoxin Level:** <1.0 EU per 1µg (determined by the LAL method).

**Applications:** SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted Molecular Mass:** 25.9kDa

**Predicted isoelectric point:** 7.1

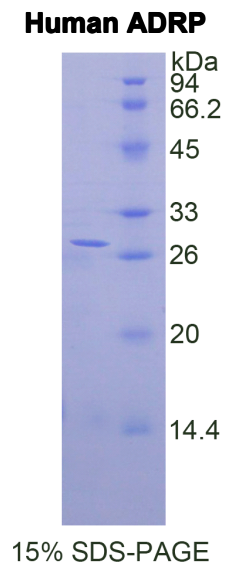



### **[ PREPARATION ]**

Reconstitute in sterile PBS, pH7.2-pH7.4.

### **[ STORAGE AND STABILITY ]**

**Storage: Avoid repeated freeze/thaw cycles.**

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

### **[ SEQUENCES ]**

The target protein is fused with N-terminal His-Tag, its sequence is listed below.

MGHHHHHSGSEF- AVDPQP SVVTRVVNLP LVSSTYDLMS SAYLSTKDQY PYLKSVCEMA  
ENGVKTITSV AMTSALPIIQ KLEPQIAVAN TYACKGLDRI EERLPILNQP STQIVANAKG  
AVTGAKDAVT TTVTGAKDSV ASTITGVMDK TKGAVTGSVE KTKSVVSGSI NTVLGSRRMMQ  
LVSSGVENAL TKSELLVEQY LPLTEEELEK EAKKVEGF DL VQKPSYYVRL GSLSTKLHSR  
AY

Gentaur Molecular Products  
Voortstraat 49  
1910 Kampenhout, BELGIUM  
Tel 0032 16 58 90 45 | Fax 0032 16 50 90 45  
[www.gentaur-worldwide.com](http://www.gentaur-worldwide.com)  
[info@gentaur.com](mailto:info@gentaur.com)