



**rP90797Ra01 50µg**  
**Angiotensinogen (AGT)**  
**Organism: Rattus norvegicus (Rat)**  
***Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

7th Edition (Revised in May, 2013)

### [ **PROPERTIES** ]

**Residues:** Asp25~Val477 (Accession # P01015),  
with two N-terminal Tags, His-tag and S-tag.

**Host:** *E. coli*

**Subcellular Location:** Secreted.

**Purity:** >95%

**Endotoxin Level:** <1.0EU per 1µg  
(determined by the LAL method).

**Formulation:** Supplied as lyophilized form in PBS,  
pH7.4, containing 5% sucrose, 0.01% sarcosyl.

**Predicted isoelectric point:** 5.4


**Predicted Molecular Mass:** 55.3kDa

**Applications:** SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

### [ **USAGE** ]

Reconstitute in sterile PBS, pH7.2-pH7.4.

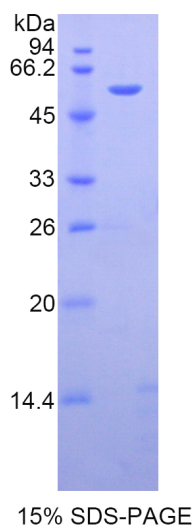

## [ STORAGE AND STABILITY ]

### **Storage: Avoid repeated freeze/thaw cycles.**

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

## [ SEQUENCES ]

The target protein is fused with two N-terminal Tags, His-tag and S-tag, its sequence is listed below.

MHHHHHSSG LVPRGSGMKE TAAAKFERQH MDSPDLGTDD DDKAMADIGS EF- DRVYIH  
PFHLLYYSKS TCAQLENPSV ETLPEPTFEP VPIQAKTSPV DEKTLRDKLV LATEKLEAED  
RQRAAQVAMI ANFMGFRMYK MLSEARGVAS GAVLSPPALF GTLVSFYLGSLDPTASQLQV  
LLGVPVKEGD CTSRLDGHVK LTALQAVQGL LVTQGGSSSQ TPLLQSTVVG LFTAPGLRLK  
QPFVESLGPFP TPAIFRSLD LSTDPVLAQK KINRFVQAVT GWKMNPLEG VSTDSTLFFN  
TYVHFQGKMR GFSQLTGLHE FWVDNSTS SVS VPMLSGTG NF QHWSDAQNNF  
SVTRVPLGES VTLIIQPQC ASDLDRVEVL VFQHDFTWI KNPPPRAIRL TLPQLEIRGS  
YNLQDLLAQA KLSTLLGAEA NLGKMGD TNP RVGEVLNSIL LELQAGEEEQ PTESAQQPGS  
PEVLDVTLSS PFLFAIYERD SGALHFLGRV DNPQNVV

## [ REFERENCES ]

1. Alvin Z., *et al.* (2011) Can. J. Physiol. Pharmacol. 89:206-215.
2. Shen Y., *et al.* (2011) Zhong Xi Yi Jie He Xue Bao 9:313-319.
3. Gard P.R., *et al.* (2011) J. Pharm. Pharmacol. 63:565-571.
4. Stobdan T., *et al.* (2011) J Renin Angiotensin Aldosterone Syst 12:93-101.



Gentaur Molecular Products  
Voortstraat 49  
1910 Kampenhout, BELGIUM  
Tel 0032 16 58 90 45 | Fax 0032 16 50 90 45  
[www.gentaur-worldwide.com](http://www.gentaur-worldwide.com)  
[info@gentaur.com](mailto:info@gentaur.com)